Walcora Rollator

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INSTRUCTION MANUAL

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Deutsch S. 1-9

English P. 10-18

Polska S. 19-27

Tedesco pp. 28-36

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EN: For colors: black, blue, red

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English P. 10-18

Table of Contents:

I. Product Details	11
II. Safety Guidelines	11
III. Specifications.....	12
IV. Rollator Components List	12
V. Assembly Instructions	14
VI. Operation Instructions	16
VII. Maintenance	18
VIII. After Sales Service	18

Thank you for purchasing a Walcora Rollator!

The Walcora Rollator, made from high-quality aluminum, is designed with safety and reliability in mind. Enhancing mobility and independence, it reflects our commitment to helping you

Focus on Life.

EN

Please carefully review these instructions before using the product and keep them for future reference. For any questions or suggestions, feel free to contact us at: service@walcora.com

Product: Rollator

Models: WA-RL-R80-RED, WA-RL-R80-BLU, WA-RL-R80-BLK, TRA38

SUNGO Europe B.V., Fascinatio Boulevard
522, Unit 1.7, 2909VA Capelle aan den
IJssel, Netherlands, Tel: +31 (0) 10 3034500

Manuf.: Foshan Oscar Medical Instrument Co. Ltd, No.2,
(Workshop C), Nanhai National Eco-Industrial
Demonstration Park, Danzao Town, Nanhai
District, Foshan City, Guangdong Province, China

Importer: R.E.S. GmbH
Röntgenstr. 15, 66482 Zweibrücken, Germany

Should a serious incident occur, we kindly ask you to report it immediately to the manufacturer as well as to the competent authority of the country you are in. Information about the competent authority can usually be found on the official websites of your country's health or medical device authorities.

Keep plastic packaging materials away from children to avoid suffocation hazards.

I. PRODUCT DETAILS

● **Intended purpose:**

The Walcora Rollator is a medical device designed to support individuals with limited mobility. The rollator improves stability and safety while walking and promotes independent movement both indoors and outdoors.

● **Intended usage:**

The Walcora Rollator is a mobility aid designed for walking and resting, suitable for indoor and outdoor use on flat surfaces. It has been developed for individuals with limited mobility who are capable of safely operating the rollator. Users must read this manual carefully to avoid injuries and minimize damage to the rollator or the surroundings. If you have any uncertainties, please contact your dealer or physician.

● **Important Note:**

The guarantee is valid only if the product is used as intended and according to the specified conditions.

● **Contraindication:**

The Walcora Rollator is not suitable for individuals with very weak arm strength, significant cognitive impairments, or poor balance. Usage should be assessed and guided by a doctor or therapist.

II. SAFETY GUIDELINES

● **Before use:**

- Always keep your fingers away from the folding scissor mechanism to avoid injuries during folding or unfolding.
- Fully unfold and lock the rollator into position before use.
- Ensure the rollator is fully opened, with the folding scissor arm straight, before use.
- Inspect the rollator regularly to ensure all screws and fasteners are secure. Tighten them if necessary.
- If you are unsure about adjusting the handle height, seek guidance from a healthcare professional or qualified dealer.
- Adjust both push handles to the same height for comfortable and safe use.
- Check that the brakes and backrest are in proper working condition before every use.
- Only use the rollator on flat, stable surfaces. Avoid uneven or sloped ground to prevent accidents.

● **Accident Prevention:**

- Do not use the rollator if it is damaged, broken, or missing any parts.
- Avoid leaning on the backrest while sitting, as it is not designed to support full body weight.
- Always use both brakes together to stop safely. Do not sit on the rollator unless the parking brake is engaged and the rollator is on a flat surface.
- Do not attempt to replace wheels or castors on your own. Contact an authorized repairer if replacements are needed.
- Maintain a safe distance from other people or obstacles to prevent collisions.
- Be cautious when crossing thresholds, railway tracks, or uneven surfaces, as the wheels can get stuck and cause a sudden stop.
- Do not use the rollator near flammable or combustible materials, or open flames.
- Keep all four wheels in contact with the ground to maintain stability and prevent accidents.

● **Restrictions of Use**

- Never stand on the rollator.
- The net bag is designed to carry up to 11 LBS (5 kg). Do not overload it.
- Do not exceed the maximum user weight of 300 lbs (136 kg).

III. SPECIFICATIONS

Net Weight
5,95 kg*

Overall length
675 mm

Overall width
640 mm

Seat size
479 mm x 240 mm

Weight capacity
136 kg

Seat height
546 mm

Handle height
890 mm – 995 mm

Packing dim
440x195x605 mm

IV. ROLLATOR COMPONENTS LIST (1)

*Without bag and accessories

IV. ROLLATOR COMPONENTS LIST (2)

Inspection on Delivery

Upon delivery of your rollator, carefully inspect the contents to ensure all items are included and undamaged. If any parts are missing or damaged, please contact our customer service team immediately at service@walcora.com

Front wheels × 2

Cane holder × 1

Main rollator frame × 1

Tray × 1 (optional)

Storage bag × 1


Upper part with brake system and rear wheels × 1

SCAN
TO VIEW THE INSTALLATION
VIDEO

V. ASSEMBLY INSTRUCTIONS

1. FRONT WHEELS

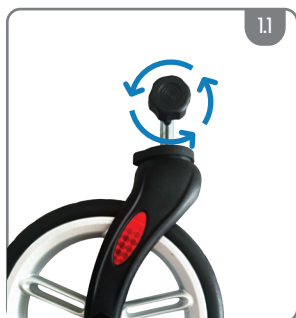1.1 Loosen the pre-assembled handwheel screws.

1.2 Insert the front wheel into the main tube of the frame, align the holes and insert the handwheel screw.


1.3 Make sure the handwheel screw is correctly seated in the hexagonal slot and tighten it securely.


2. REAR WHEELS

2.1 Loosen the pre-assembled handwheel screws.

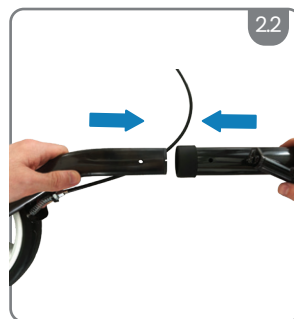2.2 Insert the rear wheel tube into the main tube of the frame.

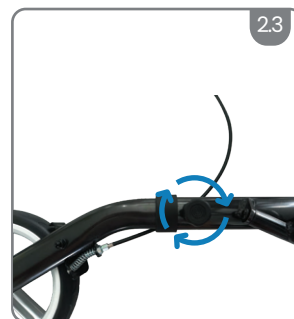2.3 Align the holes, insert the handwheel screw and tighten it securely.

The brake mechanism of the rear wheels must be on the outside of the rollator. Make sure that the handwheel screw is on the inside of the frame.

EN

V. ASSEMBLY INSTRUCTIONS

3. UNFOLDING THE MAIN FRAME

3.1 Stand the rollator upright. Then push both sides of the seat down until the frame is completely unfolded.

Attention: The rollator is only fully unfolded when a distinct “click” sound is heard, confirming the locking mechanism.

4. ATTACHING THE HANDGRIPS

4.1 Insert the handle tube into the main frame tube.

4.2 Press and hold the push button on the handle and push the handle in until the button audibly clicks into the first hole.

1. Make sure the thumb rest on the handle is installed on the inner side of the rollator.
2. Check that the button is fully engaged.

5. ATTACHING THE BAG

5.1 Hang the storage bag onto the holders of the seat frame. If the bag is not being used, it can be removed.

V. ASSEMBLY INSTRUCTIONS

6. INSTALLING THE CANE HOLDER

6. Position the cane holder over the designated slot and slide it down until it is securely fixed.

VI. OPERATING INSTRUCTIONS

1. SITTING


- The rollator provides support while walking. For example, if you want to take a break, you can sit on it.
- Place the rollator on a firm, flat surface. Activate the brake in the parking position. See section 4 (Brake system).
- The rollator must stand securely and must not slip. Sit upright.

EN

2. HANDLE HEIGHT ADJUSTMENT

If the handles are not at a comfortable height, you can easily adjust them:


- Stand upright and let your arms hang loosely: the handles should be at the height of your wrists.
- To change the height, press and hold the push button. Move the handle to the desired position and release the button so that it locks into the nearest hole (see Figure 2.1).
- Finally, check that the button is fully locked in place before using the rollator.

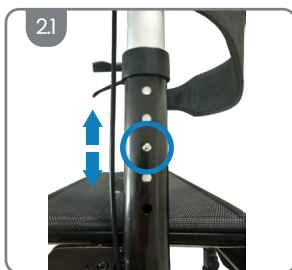3. BRAKE ADJUSTMENT

If the brakes are too loose or too tight, you can adjust them as follows:

1. The adjustment nut for the brake is located at the end of the brake cable, directly above the rear wheels on both sides.
2. First, loosen the locking nut with a 10mm spanner by turning it counterclockwise.
3. Use an 8mm spanner to adjust the brake tension. Turn the adjustment nut counterclockwise to tighten the brake cable and increase the braking force.
4. After adjustment, check the brake function carefully.

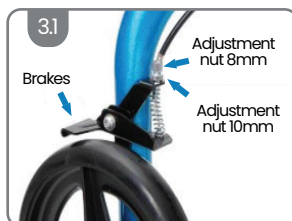VI. OPERATING INSTRUCTIONS

4. BRAKE SYSTEM

To slow down: Pull the brake lever upwards to reduce speed.

To park: Push the brake lever down until you hear a click – this locks the rear wheels. To release the brakes, pull the lever back up.


Using the seat: Always ensure the brakes are in parking position before sitting down. Test this by gently pushing the rollator to ensure it stands securely.

Notes:


- Improperly adjusted brakes can lead to injuries.
- When the brake lever is released, there must be a small gap between the brake mechanism and the wheel. See Section 3 – Brake Adjustment.

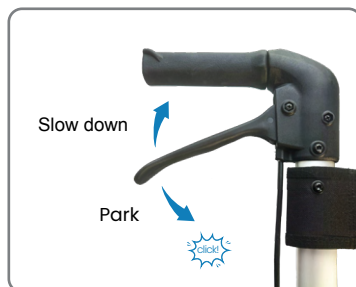5. FOLDING

5.1 Grab the handle in the middle of the seat and pull it upwards. This automatically folds the seat upward.

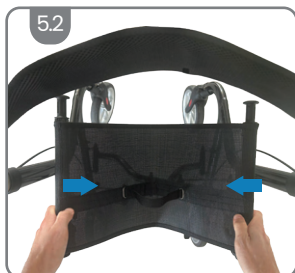5.2 Hold the rollator firmly and press both sides of the seat inward until the rollator is fully folded.

5.3 Secure the folded rollator by snapping the locking clip at the seat into place. This keeps the rollator closed and makes it easy to transport or store.

VII. MAINTENANCE

Inspection

Over time, brake cables may stretch. If adjustments are no longer possible, the cables should be replaced. Consult a qualified technician or repair service for assistance with replacing the brake cables. The castors should be checked to ensure they are secure each time before using the rollator. Additionally, check that the brakes operate correctly and that the handles and hand knobs are secure. If any part of the rollator appears loose or insecure, stop using it until it has been inspected and repaired by a qualified professional. The following shows the maintenance schedule:

Weekly inspection

- Inspect the brakes for malfunctions or wear.
- Check wheels and castors for loose bolts, nuts, or knobs.
- Ensure the frame is stable and free from cracks or visible damage.
- Verify that the handles are securely tightened and properly adjusted.
- Test the locking mechanism to ensure it functions smoothly.
- Clean any dirt or debris from moving parts, including wheels and brakes, to prevent wear.

Periodical inspection

- Inspect hand grips for wear or looseness.
- Ensure all frame bolts are secure and properly tightened.
- Inspect the front and rear axle bolts to confirm they are secure.
- Check tires for wear, cracks, or other signs of damage.
- Verify the brake cables for fraying or excessive stretching.
- Verify the folding mechanism operates smoothly and locks firmly.
- Inspect the seat and backrest (if applicable) for wear, tears, or looseness.
- Ensure accessories like baskets or cane holders are securely attached.

EN

Cleaning

Use a damp cloth or sponge with warm soapy water to clean your rollator. Domestic anti-bacterial cleaners may also be used if necessary. Make sure to remove any cleaning residue from the hand grips and seat (if applicable) to avoid them becoming slippery during use. Avoid using scouring pads or scrapers, as these can damage the surface coating.

VIII. AFTER SALES SERVICE

Replacement of Parts

Wearable parts like tires may need replacement due to normal use. Contact our customer service team for assistance with replacements.

Warranty

Your rollator includes a 2-years limited warranty covering manufacturing defects in materials and workmanship from the date of purchase. Normal wear and tear, accidental damage, and issues caused by improper use or unauthorized modifications are not covered. Warranty terms comply with local consumer protection laws.

How to Get Service

For service, contact our customer support team. Proof of purchase (e.g., a receipt) may be required for warranty claims or service requests.

Contact

For service and support, email us at: service@walcora.com

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Producent: Foshan Oscar Medical Instrument Co. Ltd, No.2,
(Workshop C), Nanhai National Eco-Industrial Demon-
stration Park, Danzao Town, Nanhai District, Foshan City,
Guangdong Province, China

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